

MedlineExamGloves The best choice for quality protection



Choosing the right glove is a challenge; today you have more choices than ever before, thanks to advanced glove technology. Medline wants to make it simple for you to evaluate those options. This new guide to Medline exam gloves clarifies the differences between major glove types to help you make an informed decision. It even contains a "Gloves by Department" section that matches glove properties to the typical requirements of specific areas in the hospital.

Medline is focusing on the future by developing new gloves that use the very latest technology to provide outstanding fit, comfort, and performance. Be sure to read about our Aloetouch® gloves, which have an aloe vera coating that soothes healthcare workers' hands, and our improved nitrile gloves, which offer great strength and are more comfortable than ever.

As always, we welcome your suggestions.

Andy Mills

President

Medline Industries, Inc.

CONTENTS

Aloetouch Gloves	3
Generation Pink Gloves	5
Nitrile Gloves	7
Premium Vinyl Gloves	9
Gloves by Department	11
Glove Properties	13
Item Index	14

Aloetouch® Gloves

The only exam gloves coated with aloe vera to help soothe skin

Every day, millions of healthcare workers trust Aloe*touch* exam gloves to protect their hands like no other glove can. Coated on the inside with pure aloe vera gel — a natural emollient — Aloe*touch* gloves soothe the skin and help to keep it soft.



Benefits of Better Skin Condition

- Greater compliance with handwashing and glove protocols
- Less chance of cross-contamination
- Higher job satisfaction among clinicians
- Potential savings from reduced turnover

Aloetouch Options

Aloetouch Ice powder-free nitrile

Coated on the inside with pure organic aloe to help soothe the skin, this glove offers the protection of nitrile and a latex-like feel all in one exam glove.

Aloetouch Ultra IC powder-free stretch vinyl

Soothes skin with natural aloe; strong and durable stretch synthetic option. Enhanced inner coating for easier donning.

Aloetouch 3G powder-free stretch vinyl

Coated with skin-soothing aloe — the best-performing vinyl glove on the market.

Aloetouch powder-free nitrile

Promotes and maintains soft skin — the industry's best choice for complete protection. Available in both 9.5-inch and 12-inch lengths.

Aloetouch powder-free latex

Soothes skin, with low protein levels to help reduce latex sensitivity.

Thanks to Medline for providing a glove with an aloe vera coating that truly soothes the skin. As an infection control practitioner who promotes excellent skin care for healthcare workers, I truly believe the aloe glove enhances hand care.

Louise Standridge, RN Director of Epidemiology, Carraway Methodist Medical Center

Visit www.aloetouch.com to learn more!











aboratory







Generation Pink® Gloves

One Mission. Three Gloves.

Medline started the Generation Pink® movement with a single vision: A World Without Breast Cancer is in Our Hands.TM It's now the largest philanthropic cause in our company, where we partner with the National Breast Cancer Foundation to support free mammograms for underserved women, build awareness at forums for nurses and other healthcare workers, and even spread awareness through social media and events such as The Pink Glove DanceTM (PinkGloveDance.com).



Medline now has three Generation Pink gloves that help financially support programs for breast cancer awareness, early detection and prevention among underserved women. Wearing the pink gloves helps show your support.

Pink Glove Options

Generation Pink Pearl® powder-free nitrile

Part of the Aloetouch collection, Pearl gloves have soothing aloe on the inside and an attractive shimmer on the outside. Textured throughout. For each case you purchase, \$1 is donated to the National Breast Cancer Foundation.

Generation Pink® powder-free nitrile

Excellent tactile sensitivity. Textured fingertips. For each case you purchase, \$2 is donated to the National Breast Cancer Foundation.

Generation Pink® powder-free 3G vinyl

The most advanced vinyl available. For each case you purchase, \$1 is donated to the National Breast Cancer Foundation.

















Nitrile Gloves

Trust our nitrile gloves for strength with a soft touch

Medline has invested heavily in improving the technology of nitrile gloves to increase the softness while maintaining excellent chemical resistance and barrier protection. Our nitrile gloves are now more comfortable to wear, and new manufacturing methods have made them more affordable.





Proven Chemical Resistance

Nitrile gloves are stronger and possess greater chemical resistance than both latex and vinyl gloves. They are generally recommended for use in areas that may involve higher levels of stress on exam gloves (e.g., ER and trauma units) or areas that deal with harsh chemicals or liquid immersion (e.g., Lab, Oncology, Environmental Services, Endoscopy, and Pharmacy).

Nitrile Options

SensiCare Ice® powder-free nitrile

Provides great sensitivity and dependable protection. Chemo-tested.

Curad® powder-free nitrile

Excellent strength and high chemical resistance for hospital-wide protection. Chemo-tested.

SensiCare® powder-free nitrile

High-performance latex-free protection. A soft, strong glove with excellent tactile sensitivity. Chemo-tested.

Accutouch® powder-free nitrile

Comfortably conform to the hand shape. Textured throughout. Tight cuffs for extra protection. Chemo-tested.

SensiCare Silk® powder-free nitrile

Our highest level of tactile sensitivity. Unique dark-blue color.

SensiCare Free™ powder-free nitrile

Low potential for dry, irritated skin. Made without the accelerants that are often used when producing nitrile and latex. Chemo-tested.

MediGuard® powder-free nitrile

Extra-sensitive nitrile glove; economical choice. Chemo-tested.

SensiCare 12" powder-free nitrile

Extended cuffs for extra protection. Excellent tactile sensitivity. Textured fingertips. Chemo-tested.

For Aloetouch nitrile options, see page 4.
For Generation Pink nitrile options, see page 6.















Premium Vinyl Gloves

Medline vinyl synthetic gloves: all-around performance made better by technology

Advances in vinyl synthetic glove technology have enabled Medline to make significant improvements in the features that matter most to you: strength, elongation, tactile sensitivity, and comfort. Our third-generation vinyl synthetic gloves are nothing like the vinyl synthetic gloves of a few years ago, thanks to advanced-formulation polymers.



Ideal for Latex-Free Environments


Today's improved vinyl synthetic gloves are a cost-effective, latex-free alternative well suited to the majority of clinical applications. Except for a few instances that call for specialty gloves (e.g., chemotherapy gloves in Oncology), stretch vinyl gloves are an excellent choice for economical latex-free barrier protection.

Three Generations... What's the Difference?

- First: traditional vinyl material containing PVC resin and plasticizers.
- Second: stretch vinyl material with softer PVC resins for better fit, feel, and performance.
- Third: stretch vinyl material containing ultra-highmolecular-weight resins and special softening and strengthening agents; offers strength, comfort, and barrier protection other vinyl gloves cannot.

Comparing performance:

These charts show how tensile strength and elongation improve with each successive generation.

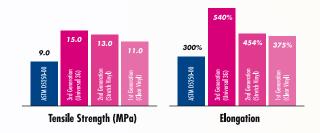

Vinyl Synthetic Options

Curad 3G third-generation powder-free stretch vinyl

Outstanding comfort and performance. Proprietary formulation offers remarkable softness, strength, and stretch ability.

Curad second-generation powder-free stretch vinyl Improved material offers a good fit and tactile sensitivity.

For Aloetouch vinyl options, see page 4. For Generation Pink vinyl options, see page 6.











Gloves by Department



Critical and Noncritical Care Units

Units such as Intensive Care, Med-Surg, Radiology, Respiratory, and Physical Therapy perform procedures that do not generally place significant levels of stress on examination gloves. Facilities have greater flexibility in selecting an exam glove for these areas. Powder-free nitrile, powder-free latex, and second- or third-generation stretch vinyl gloves are all excellent choices.

Special Requirements	Recommended Glove Type			
Limited	Powder-free nitrile Powder-free latex 2nd- or 3rd-generation stretch vinyl			



Emergency Room

Exam gloves selected for use in the Emergency Department should provide protection against bloodborne pathogens and be strong enough to withstand the rigors of the tasks performed. Nitrile and latex gloves are the recommended choices for high-stress activities. High-performance third-generation vinyl gloves can be used as a complement for routine patient care.

Special Requirements	Recommended Glove Type
Durable gloves that can be used for all procedures	Powder-free nitrile or powder-free latex



(a) Labor and Delivery

Labor and Delivery generally prefers to use latex-free gloves to avoid any potential latex allergy issues with newborns. Sterile gloves are encouraged for actual deliveries, which may require more invasive procedures. Medline offers both stretch vinyl and nitrile sterile glove options.

Special Requirements	Recommended Glove Type			
Latex-free	Powder-free nitrile or 3rd-generation stretch vinyl			
Sterile procedures	Sterile powder-free nitrile or sterile 2nd-generation stretch vinyl			



Laboratory

Laboratory clinicians have a variety of glove needs. Some technicians deal with harsh chemicals and require higher chemical resistance, while others have greater grip requirements and need a glove that will perform well in both wet and dry situations. All clinicians require comfortable, high-performing gloves. Those who deal with harsh chemicals should use nitrile gloves, but latex and third-generation high-performing vinyl gloves are appropriate for all other functions.

Special Requirements	Recommended Glove Type			
Strong chemical resistance	Powder-free nitrile			
Greater grip requirements Comfort Tactile sensitivity Durability	Powder-free latex 3rd-generation stretch vinyl			



Phlebotomy

Phlebotomists need a glove with excellent tactile sensitivity that allows them to find a vein hidden beneath the skin.

Special Requirements	Recommended Glove Type
Excellent tactile sensitivity	2nd- or 3rd-generation stretch synthetic Powder-free latex Powder-free nitrile



Pharmacy/Oncology

Clinicians who are preparing, compounding, handling, or administering harsh chemotherapy drugs in these departments should use chemotherapy-tested nitrile gloves.

Special Requirements	Recommended Glove Type			
Resistance to chemotherapy agents	Chemotherapy-tested powder-free nitrile			

Environmental Services

Environmental Services employees may handle harsh chemical cleaning compounds and ideally should use nitrile gloves, which provide the most effective chemical resistance of all glove materials.

Special Requirements	Recommended Glove Type			
Strong chemical resistance	Powder-free nitrile			



Glove Properties

Viral Penetration

The viral penetration test determines if a glove provides an effective barrier against one of the smallest measurable viruses in the world, phi-X174. If the glove is a barrier to the phi-X174 virus, it also provides a barrier against larger viruses.

All of Medline's exam gloves have passed the viral penetration test.

Chemical Resistance and Barrier Guide

	-	Curad Stretch Vinyl		d 3G ı Vinyl		rad Nitrile	Aloe <i>touch</i> 12" Chemo Nitrile	
	Permeation Rate	Breakthrough (min.)	Permeation Rate	Breakthrough (min.)	Permeation Rate	Breakthrough (min.)	Permeation Rate	Breakthrough (min.)
Acetaldehyde	7.31	16	6.84	18	3.43		1.85	240
Ammonium Hydroxide 30-70%	1.19	335	1.26	349	2	391	0.74	410
Carmustine	N/A	N/A	N/A	N/A	0.1	480	0.1	480
Citric Acid 30-70%	0.26	358	0.1	366	0.1	367	0.1	394
Chlorox Bleach 10%	16.8				6.1	219	4.4	286
Cyclohexanol	14.7		13.6		N/T	N/T	6.2	382
Ethanol	17.1	5	16.3	6	8.6	146	4.3	216
Formaldehyde 30-70%			3.6		2.81	351	1.64	388
Formic Acid	9.2	365	7.3	378	16.2	247	10.4	286
Freon 113 or TF	23.1	12	24.6	13	14.8	360	8.6	380
Glutaraldehyde	N/T	N/T	0.18	480	0.1	480	0.1	480
Hydrazine	10.8	352	13.6	364	21.6	248	11.4	376
Hydrogen Peroxide 30%	15.1	253	16.3	248	5.8	360	4.9	360
Isopropyl Alcohol 70%	12.1		13.8	128	4.3	378	3.8	377
Isopropyl Alcohol 90%	7.1	150	6.2	156	4.3	378	2.14	389
Methanol	52.3	4	43.8	4.6	29		21	124
Potassium Hydroxide	20.2	235	18.4	246	13.2	273	9.3	386
Povidone Iodine 10%	7	280	6.8	283	23.7	253	16.4	275
Quaternary Cleaners A5 by 3M	N/T	N/T	0.82	53	13.8	140	0.9	173
Triclosan	0.03	266	0.03	342	N/T	N/T	N/T	N/T

Rating Code: The basis of the rating code is breakthrough time in minutes only. We have color-coded the ratings for your convenience. The following table explains the rating system.

Not Recommended
Breakthrough occurs in
less than 30 minutes.
Avoid use of the glove
with this chemical.

Caution
Breakthrough occurs
in 30 minutes to
2 hours. This glove
is suitable for this
chemical under careful
control of its use.


Recommended
Breakthrough occurs
after 2 hours. This
glove is well suited
for application with
this chemical.



Permeation rate measured in $\mu g/cm^2/minute$.

ASTM F739, Standard Method for Determining Resistance to Chemical Permeation Under Conditions of Continuous Contact. Testing completed by Lowenkamp R&D Laboratories Inc. Test reports on file.

Tensile Strength and Elongation

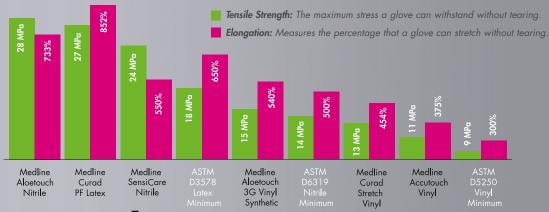

Item Index

	Description	Quantity	x-small	small	medium	large	x-large
	Aloe <i>touch</i> Ice PF Nitrile	2000/case	MDS195283	MDS195284	MDS195285	MDS195286	MDS195287
	Aloe <i>touch</i> Ultra IC PF Vinyl	1000/case	MDS195073	MDS195074	MDS195075	MDS195076	MDS195077
	Aloetouch 3G PF Vinyl	1000/case	MDS195173	MDS195174	MDS195175	MDS195176	MDS195177
	Aloe <i>touch</i> PF Nitrile	1000/case	MDS195083	MDS195084	MDS195085	MDS195086	MDS195087
Aloe	Aloetouch 12" PF Nitrile	500/case	_	MDS195184	MDS195185	MDS195186	MDS195187
	Aloe <i>touch</i> PF Latex	1000/case	CUR8153	CUR8154	CUR8155	CUR8156	CUR8157
	Aloetouch 12" PF Nitrile Singles, Sterile	400/case	_	MDS194081	MDS194082	MDS194083	MDS194084
	Aloe <i>touch</i> 12" PF Nitrile Pairs, Sterile	200pr/case	_	MDS194085	MDS194086	MDS194087	MDS194088
	Aloetouch 9" PF Nitrile Singles, Sterile	400/case	_	_	MDS194285	MDS194286	_
	Generation Pink Pearl PF Nitrile	1000/case	PINK5083	PINK5084	PINK5085	PINK5086	PINK5087
Pink	Generation Pink PF Nitrile	2000/case	PINK6800	PINK6801	PINK6802	PINK6803	PINK6804
	Generation Pink PF Vinyl	1000/case	PINK6073	PINK6074	PINK6075	PINK6076	PINK6077
	SensiCare Ice PF Nitrile	2000/case	486800	486801	484802	486803	486804
	Curad PF Nitrile	1000/case	_	CUR8314	CUR8315	CUR8316	CUR8317
	SensiCare PF Nitrile	1000/case	484800	484801	484802	484803	484804
a)	Accutouch PF Nitrile	1000/case	_	MDS192084	MDS192085	MDS192086	MDS192087
Nitrile	SensiCare Silk PF Nitrile	2000/case	MDS7083	MDS7084	MDS7085	MDS7086	MDS7087
	SensiCare Free PF Nitrile	2000/case	MDS19083	MDS7084	MDS19085	MDS19086	MDS19087
	MediGuard PF Nitrile	2000/case	MG6050	MG6051	MG6052	MG6053	MG6054
	SensiCare 12" PF Nitrile	500/case	_	MDS1284	MDS1285	MDS1286	MDS1287
	PF Nitrile Pairs, Sterile	200pr/case	_	MDS198314	MDS198315	MDS198316	MDS198317
	Curad 3G PF Vinyl	1000/case	CUR8233	CUR8234	CUR8235	CUR8236	CUR8237
	Curad PF Stretch Vinyl	1000/case	_	CUR8224	CUR8225	CUR8226	CUR8227
	Ultra PF Stretch Vinyl	1000/case	_	MDS193074	MDS193075	MDS193076	MDS193077
Viny	Ultra 12" PF Stretch Vinyl	1000/case	_	MDS191274	MDS191275	MDS191276	MDS191277
	Accutouch PF Vinyl	1000/case	MDS192073	MDS192074	MDS192075	MDS192076	MDS192077
	SensiCare PF Vinyl Singles, Sterile	400/case	_	484401	484402	484403	_
	SensiCare PF Vinyl Pairs, Sterile	200pr/case	_	484405	484406	484407	-
	Curad PF Latex	1000/case	CUR8103	CUR8104	CUR8105	CUR8106	CUR8107
	MediGuard PF Latex	1000/case	_	MG8701	MG8702	MG8703	MG8704
ex	MediGuard PWD Latex	1000/case	MG1203	MG1204	MG1205	MG1206	MG1207P
Latex	Eudermic® MP PF Latex	500/case	_	485601	485602	485603	485604
	PF Latex Singles, Sterile	400/case	_	MDS194134	MDS194135	MDS194136	_
	PF Latex Pairs, Sterile	200pr/case	_	MDS194034	MDS194035	MDS194036	_ /



Medline Industries, Inc. One Medline Place Mundelein, IL 60060

Medline United States 1-800-MEDLINE (633-5463) www.medline.com info@medline.com Medline Canada 1-800-396-6996 www.medline.ca canada@medline.com Medline México 01-800-831-0898 www.medlinemexico.com mexico@medline.com

Some products may not be available for sale in Mexico or Canada. We reserve the right to correct any errors that may occur within this brochure. ©2010 Medline Industries, Inc. Medline is a registered trademark of Medline Industries, Inc. Medline, Aloetouch, Aloetouch Ice, Curad, Eudermic, Generation Pink, MediGuard, Pearl, SensiCare, SensiCare Ice and SensiCare Silk are registered trademarks of Medline Industries, Inc. A World Without Breast Cancer is in Our Hands, Pink Glove Dance, Pink Pearl and SensiCare Free are trademarks of Medline Industries, Inc. Color Pink and Pink Nitrile are registered trademarks of Colur World, LLC.

MKT210314/LIT507/10M/TKZ22